

Patient history


Last Modified on 29/10/2024 11:57 am AEDT

A patient's History View records all information added using Zedmed's modules and clinical notes.

Each encounter starts with a **date**, **time** and **duration** in bold, and the information can be searched, filtered and actioned as explained in the sections below.

Search for specific information

Select the **Search History** button, enter a term and select **Find Next**. This Searches the range in the **Visible** field.

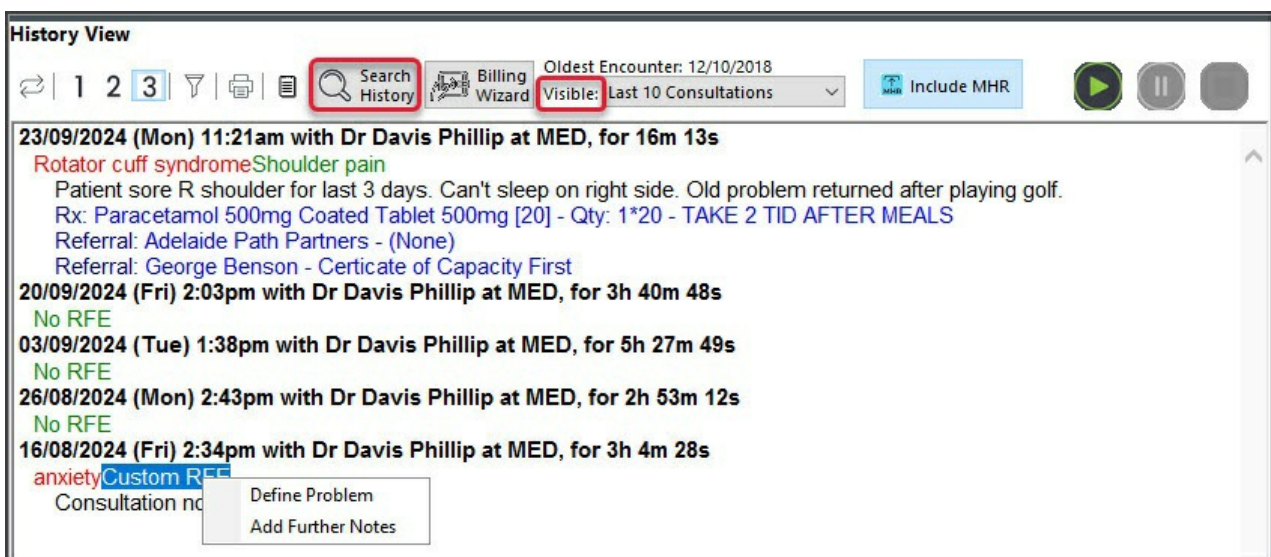

The screenshot shows the 'History View' interface. At the top, there are navigation buttons (1, 2, 3) and a search bar labeled 'Search History'. To the right, there are buttons for 'Billing Wizard' and a dropdown menu for 'Visible: Last 10 Consultations'. There is also an 'Include MHR' button and play/pause icons. The main area displays a list of encounters:

- 23/09/2024 (Mon) 11:21am with Dr Davis Phillip at MED, for 16m 13s**
 Rotator cuff syndrome **Shoulder pain**
 Patient sore R shoulder for last 3 days. Can't sleep on right side. Old problem returned after playing golf.
 Rx: Paracetamol 500mg Coated Tablet 500mg [20] - Qty: 1*20 - TAKE 2 TID AFTER MEALS
 Referral: Adelaide Path Partners - (None)
 Referral: George Benson - Certificate of Capacity First
- 20/09/2024 (Fri) 2:03pm with Dr Davis Phillip at MED, for 3h 40m 48s**
 No RFE
- 03/09/2024 (Tue) 1:38pm with Dr Davis Phillip at MED, for 5h 27m 49s**
 No RFE
- 26/08/2024 (Mon) 2:43pm with Dr Davis Phillip at MED, for 2h 53m 12s**
 No RFE
- 16/08/2024 (Fri) 2:34pm with Dr Davis Phillip at MED, for 3h 4m 28s**
 anxiety **Custom RFE**
 Consultation no []

A context menu is open over the 'Custom RFE' entry, showing options: 'Define Problem' and 'Add Further Notes'.

- Headers = **black**
- RFEs = **green**
- Problems = **red**
- Scripts & referrals = **blue**
- Notes = black

Filter what is displayed

Adjust the amount of detail shown

Select the **1**, **2** and **3** buttons to increase or decrease the level of detail shown.

Selecting **1** displays one line per appointment, and the filter must be changed to **2** or **3** to show information for that appointment.

Period shown


Select **Last 10**, **Last 1 Year**, **3 Years** or **All** in the **Visible** field's drop-down to set the period shown.

Include MHR uploads

Toggle the button on to show information added to the patient's **MyHealthRecord** by other providers.

Filter by a specific problem, medication or allergy

- a) In **Summary views**, select the **Problems, Allergies, Pregnancy** or **Medications** tab.
- b) Right-click a specific entry and select **Filter History View** to these **Encounter Segments**.

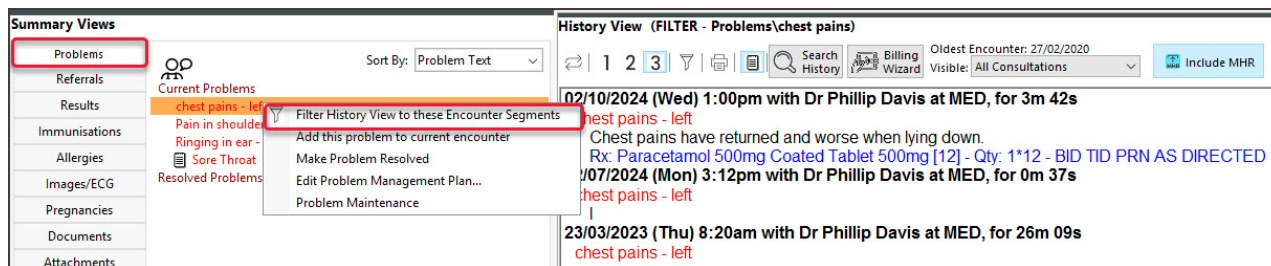

Update what is displayed

Remove notes

Up until midnight on the day the information was recorded, you can reopen the encounter. Any changes made, including the removal of clinical notes, will be reflected in the History View.

Add Notes

The day after the clinical encounter, you can append an encounter's notes by right-clicking the encounter date and selecting **Add Further Notes**. This opens a field that will add **New Addendum** notes to that encounter.

Add a Problem/RFE

To add to a previous consult, right-click the text (red or green) and select Define Problem.

Add information from the current encounter

Selecting the refresh icon will add information from the Current Encounter.

Reprint/Rename a referral

Right-click the blue referral text and select open or rename.